


## Tubular Glass Syringe systems

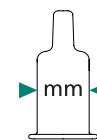

### RTF® syringe systems Luerlock

0.5 ml	47.6	6.85	4.65
1.0 ml long	54.0	8.15	6.35
1.0 ml standard	35.7	10.85	8.65
1.5 ml	43.2	10.85	8.65
2.25 ml	54.4	10.85	8.65
3.0 ml	72.2	10.85	8.65

#### Material:

Type I glass (EP/USP/JP compliant)

- RTF® – Ready-to-Fill
- System is completed by bulk or nested stoppers (see page 2.8)
- RTF® process steps: washing, siliconisation, assembly, nesting
- ETO sterilisation
- Baked-on siliconisation as option (see page 6.1)
- Drug Master File Type III
- Suited to standard filling and packaging equipment
- Customised ceramic printing with heavy-metal-free inks
- Optional: special treatment with ammonium sulphate for greater surface resistance (see page 6.2)

Tip caps in various rubber formulations and designs available on request as well as Gerresheimer TELC (Tamper-Evident Luerlock Closure)

Our technical support team would be happy to provide you with further information. Please contact [info-buende@gerresheimer.com](mailto:info-buende@gerresheimer.com).